


Clearfil® Ceramic Primer

The future of Ceramic Primer is clear

Introducing Clearfil® Ceramic Primer, the single-component advance in high-strength primers for all ceramic applications

Finally! A single component primer with the strength and durability of Kuraray's "5 Star" Porcelain Bond Activator is now available. Clearfil® Ceramic Primer is a single-component adhesive primer used to enhance the bond strength between resin-based materials and any type of ceramic restorative material – porcelain, Zirconia, Alumina, Lucite, Lithium Silicate or composites – with a simplified operation process. It contains the phosphate monomer "MDP," which bonds strongly to metal oxides and contains a Silane Coupler to bond with SiO₂-based ceramics. Clearfil® Ceramic Primer is an easy-to-use, no mix, water-free solvent that evaporates quickly. It does not form a sticky film, and a complicated HF treatment is not required. The newly discovered chemical equilibrium state between MDP and the Silane Coupler offers better shelf-life stability than other single-component products. Clearfil® Ceramic Primer provides simple and easy operation for indirect restoration, using resin cement, and for intra-oral repairs of fractured ceramic restoratives by composite.

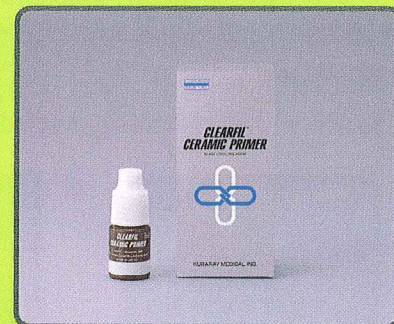


World's Leading Manufacturer of Dental Adhesive Products

**KURARAY
DENTAL**



www.kuraraydental.com

Experience the ceramic primer that is structured to be strong

Clearfil® Ceramic Primer improves handling without mixing. It's easy-to-use, evaporates quickly, and won't leave a sticky film after application.

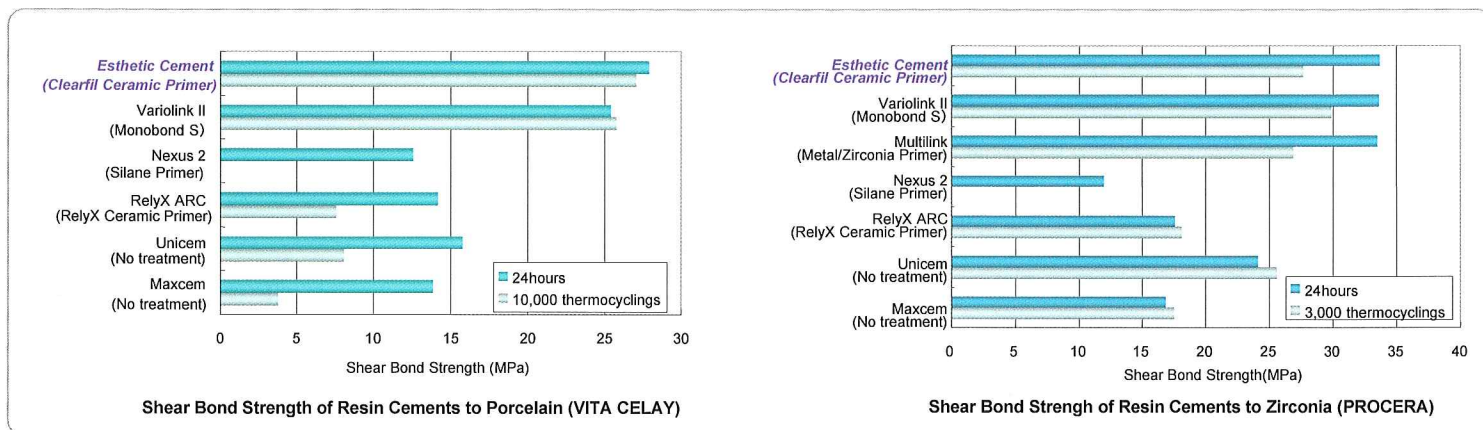


The molecular structure provides superior adhesion to various substrates with improved handling.

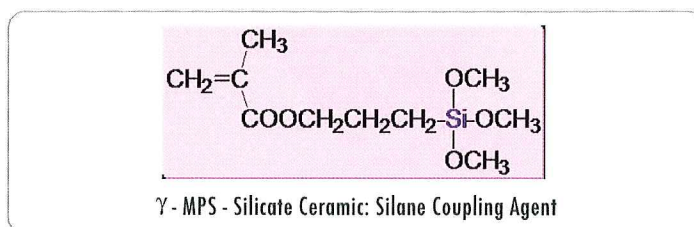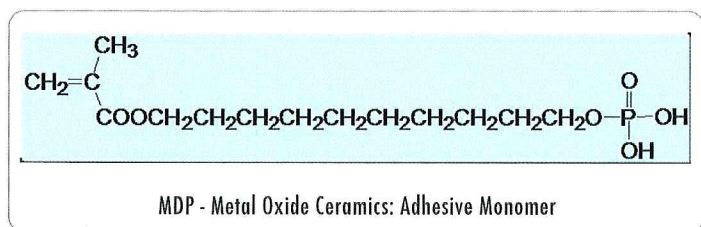

From the popular Clearfil® family come the newest developments in dental adhesives

The Clearfil® family of adhesive and restorative products is unsurpassed in reliability, esthetic beauty, and technological innovation. Sample all of the latest adhesives in the Clearfil line from Kuraray Dental, and prove it to yourself.

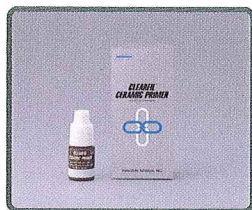

Clearfil® Ceramic Primer is a single-component adhesive primer that enhances the bond strength between resin-based materials and any type of ceramic restorative material.



Clearfil® DC Bond is a breakthrough in dual-cure adhesives that exhibits superb performance in any restorative situation, including veneers.



Clearfil® Esthetic Cement a composite-like material, has a high filler load and delivers outstanding adhesion, long wear, and high color stability.

Clearfil® Ceramic Primer Products

Item Number	Product Name	Contents
#2550KA	Clearfil® Ceramic Primer	1 Primer 4ml