



# Knit-Pak™

## Knitted Retraction Cord



A stainless steel cutting blade inside the cap automatically cuts the cord upon closing.



The cord piece is held by the cap until removed, preventing it from dropping. Protects cord piece and container from cross-contamination.

**Knit-Pak™ gingival retraction cord is a non-impregnated knitted cord made from 100% cotton.**


Knit-Pak is ideal to use in combination with Hemodent® hemostatic agent or Traxodent® for tissue displacement and management prior to taking an impression. The knitted cord contains numerous interlocking chains that carry greater quantities of hemostatic solution compared to twisted or braided cord. Knit-Pak easily packs into sulcus by compacting under pressure. It will not fray and resists popping out of the sulcus. If touched by a bur, Knit-Pak will cut clean and not get entangled.

**Knit-Pak's 6 different sizes are color coded for quick and easy identification.**

Each Knit-Pak contains 100 inches (254 cm) of cord.

### Ordering Information

REF	Size	Color	REF	Size	Color
9007551	000	Green	9007554	1	Blue
9007552	00	Brown	9007555	2	Orange
9007553	0	Purple	9007556	3	Yellow

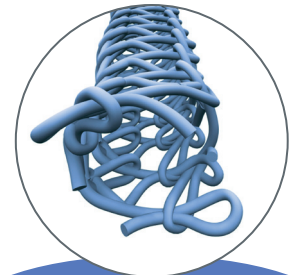

**Non-impregnated Knitted Weave**

# traxodent®

Hemodent® Paste Retraction System


**Traxodent® from Premier® provides predictable hemostasis and soft tissue management in minutes.**

## Easy, effective hemostasis and retraction


- The sleek syringe with bendable tip permits easy application.
- After two minutes it is rinsed, leaving a clean, accessible margin.

## Traxodent is gentle, absorbent and fast.

- The soft paste produces gentle pressure on the sulcus while it absorbs excess crevicular fluid.
- Aluminum chloride creates an astringent effect without irritating or discoloring surrounding tissue.

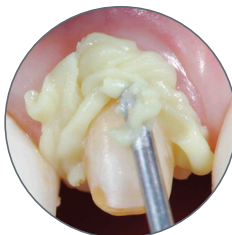

*Traxodent can be used alone.*

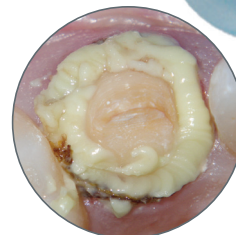

*Traxodent can replace the second cord in the double-cord technique.*



Awarded Top 100 Products in 2010 & 2011



- 9007091 **Traxodent Value Pack** (25 syringes and 50 tips)
- 9007092 **Applicator Tips** (50)
- 9007093 **Traxodent Starter Pack** (7 syringes and 15 tips)

## Premier® Retraction Caps

Cotton caps work with Traxodent to quicken retraction and hemostasis.

- Easily guides Traxodent into the sulcus when patient closes
- Compression assists and quickens hemostasis
- Anatomic form allows easy placement

9048251 **Size 1** (120)

9048253 **Size 3** (120)

9048255 **Size 5** (60)

9048260 **Assorted** - (30) #5, (15) #1, (15) #3



## Hemodent® Buffered Aluminum Chloride Epinephrine-free Hemostatic Liquid

- Effectively stops minor gingival bleeding
- Epinephrine-free to avoid cardiac reactions
- Stable and offers a long shelf-life

9007071 **10cc bottle**

9007072 **20cc bottle**

9007073 **40cc bottle**

## Premier® Cord Packers

Premier offers a selection of cord packing instruments – contact your dealer for the complete line!

- Balshi and Sherman packers, with or without serrations
- Round or taper ends in angle and straight configurations, with or without serrations



*Please contact your dealer concerning product availability in your country.*



**Premier® Dental Products Company** • 1710 Romano Drive, Plymouth Meeting, PA 19462, USA  
Phone 888-670-6100 • 610-239-6053 • Fax 610-239-6171 • www.premusa.com • E-mail: dentalinfo@premusa.com



8610325 111117 Rev2 TRI20M