


# Ultrafast Dental Laser

## Wireless is Effortless

Amazing Laser for Everyday Practice

*No More Setup Time and Messy Cords*

**3 Watts Continuous Power**

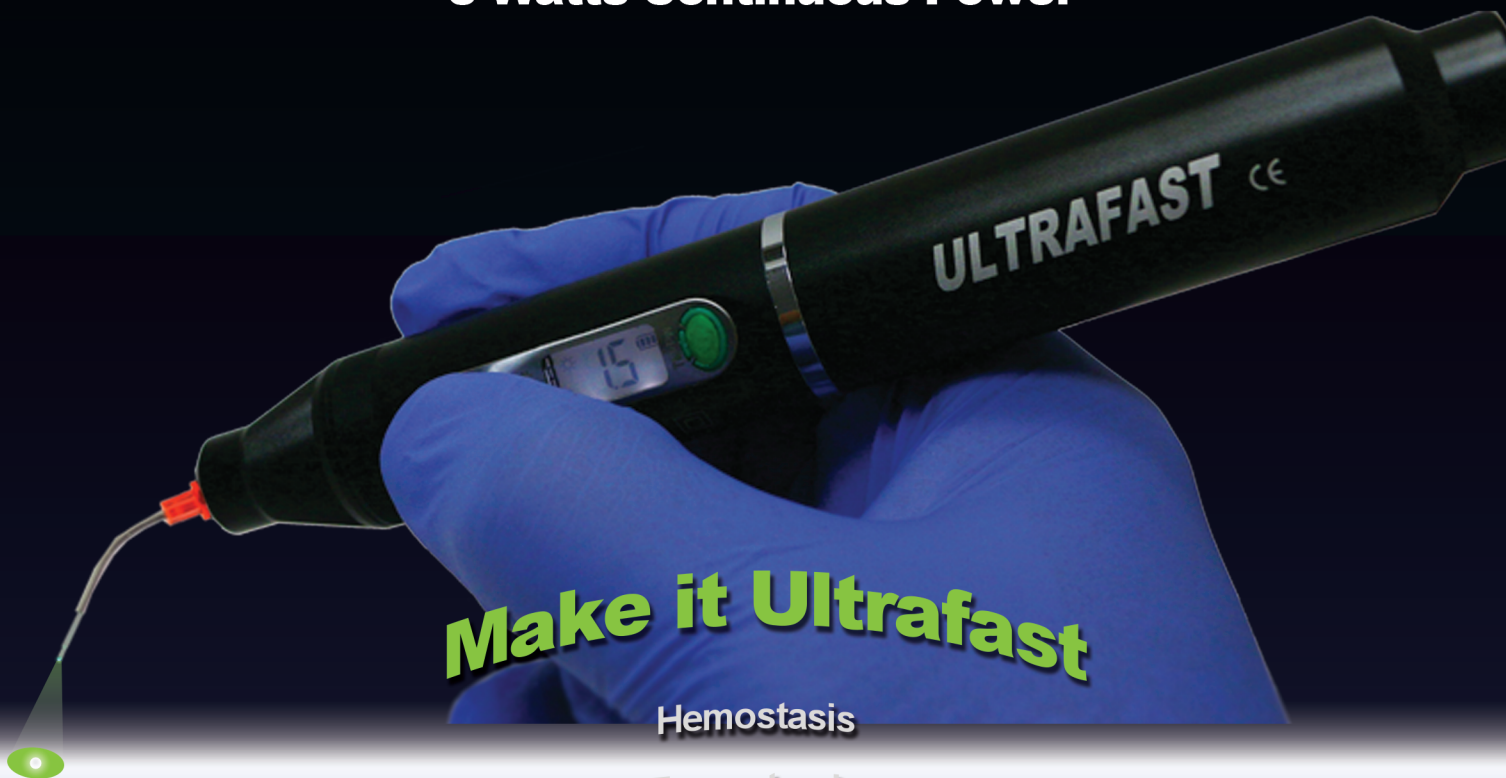

### Make it Ultrafast

Hemostasis

Frenectomies

Gingivectomies

Incision / Excision

Root Canal Therapy

Canker Sores / Aphthous ulcers

Periodontal Sulcular Debridement

Gingival Troughing for Impression

# Ultrafast Dental Laser

## Amazing Power and Style at Your Fingertip!

### Distinctive Green Aiming Beam

Easier to see treatment spot



### Sleek Wireless Design

- Patent-pending wireless laser with distinctive green aiming beam:
  - up to 3W CW Power for ultrafast soft tissue sculpting
  - up to 5W Pulsed Power for precision cut
  - Battery and thermal bars for smart control
- Eliminate cords and setup time
- 810nm wavelengths for advanced performance
- Easy slide in laser shield
- Unsurpassed 45 minutes running time
- 90 minutes fast recharge time
- Cordless with direct power plug that eliminates downtime
- Sleek metal keypad with large activation and program buttons
- Super slim and ergonomic
  - Weighs only 3.5oz
  - Measures at 6.4 inch

### Disposable Fiber Tip (self-bendable)

Save your time

Avoid cross-infection

### Built to Last

Durable scratch resistant solid metal construction distributes heat, prevent scratches, discoloration and damage from accidental falls.



### Smart Lithium-Ion Battery You Can Count On!

- Cordless with 45 minutes run time at normal cutting power
- No need to change battery in the middle of a cycle
- Multiple charge options
  - Space saving counter stand
  - Direct plug in
- Constant output throughout battery run time



Battery Assembly

Charging Stand

### PN: 7300100 Ultrafast Diode Laser Kit

Including Laser Handpiece, Charging Stand, Disposable Fiber Tips (5), Barrier Sleeves (50), Power Adapter, Laser Shield, and Laser Goggles (2).

**2-year warranty (1-year battery)**



### DENTLIGHT INC

1825 Summit Ave. #210  
Plano, TX 75074, US  
Tel: 800-763-6901  
Fax: 972-346-6550  
Email: sales@dentlight.com